

N95 Face Make

It can be used for respiratory protection of some particles, such as dust generated in the process of grinding, cleaning and handling minerals, flour and some other materials, etc., and it can also be used for liquid or non oil particles that do not produce harmful volatile gases due to spraying. It can effectively filter and purify the abnormal smell inhaled (except toxic gas), and help reduce the exposure level of some inhalable microbial particles (such as mold, anthrax, tuberculosis, etc.)

It does not eliminate the risk of contact infection, illness or death.

Product Name	N95 Face Mask
Specifications	On-ear multiplayer mask Filtration efficiency $\geq 95\%$ anti-virus / prevent flu / anti smoke / dust
Materials	Non woven, soft cotton, melt blown filter
Color	White
Executive Standard	GB19083-2010
Certificate	FDA & CE
Period of validity	Three years
Package	50 pieces in one box
Mini Order	50,000
Product Name	1. Payment : T/T 50% deposit, and the balance will be paid before delivery 2. Delivery : Within 3-10days after receiving deposit



N95
GB19083—2010



口罩结构分解图

Structure breakdown of respirator



-第一层: SS拒水无纺布白色50g/㎡
First floor: SS water repellent non-woven white 50g / m²
-第二层: 熔喷无纺布25g/㎡
Second layer: melt blown non-woven fabric 25g / m²
-第三层: 熔喷无纺布25g/㎡
The third layer: melt blown non-woven fabric 25g / m²
-第四层: SS拒水无纺布白色50g/㎡
Fourth layer: SS water repellent non-woven white 50g / m²

50 / 盒

30盒 / 箱

1500pcs/carton

重量: 15.5KG

纸箱尺寸: 62cmX56cmX63cm



CERTIFICATE

of MDD Notification

Reference number: JH-ERA-20291V35

Issued date: March 25, 2020

This certificate will be automatically void if the Notification is rejected by the EU Authorities or upon termination of the EAR.

This is certify that, According to Medical Device 93/42/EEC(MDD), we accepts the appointment to be the Authorized European Representative for products which listed in the attached agreement between below manufacturer and Luxus Lebenswelt GmbH.

Manufacturer: Henan Kanghua Pharmaceutical Co., Ltd
Address: Changdian Industrial Park, Gudun Town, Mengzhou City, Jiaozuo City, The People's Republic of China

The Manufacturer declares that the Medical device complies with de Directive including all essential requirements.

According to Medical Device 93/42/EEC(MDD), the European Databank on Medical Devices (EUDAMED) is established as of May 1, 2011, the German Competent Authority is notified of the manufacturer's Medical Devices and has allocated registration numbers shown in:


Disposable Medical Mask, UMDN code: 21-734

Protective Clothing, UMDN code: 27-725

Where the manufacturer affix the CE marking to the product listed they must ensure that all the requirements of the appropriate EU directive(s) and standards have and continue to be met.



Luxus Lebenswelt GmbH
Kochstr. 1, 47877, Willich, Germany
info.m@luxushe.de

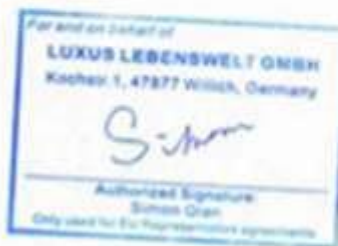



Fiscal Year 2020

CERTIFICATE OF FDA REGISTRATION

This certifies that:

HENAN KANGHUA PHARMACEUTICAL CO.,LTD
Changdian Industrial Park,Gudan Town,Mengzhou City,
mengzhou, Henan, 454761, CHINA

has completed the FDA Establishment Registration and Device Listing with the US Food & Drug Administration, through UCL-REG SERVICE INC.

Owner/Operator Number: 10065063

Listing No.	Product Code:	Device Name:
D380597	KHA	MASK,SCAVENGING

UCL-REG SERVICE INC. will confirm that such registration remains effective upon request and presentation of this certificate until the end of the calendar year stated above, unless said registration is terminated after issuance of this certificate. UCL-REG SERVICE INC. makes no other representations or warranties, nor does this certificate make any representations or warranties to any person or entity other than the named certificate holder, for whose sole benefit it is issued. This certificate does not denote endorsement or approval of the certificate-holder's device or establishment by the U.S. Food and Drug Administration. UCL-REG SERVICE INC. assumes no liability to any person or entity in connection with the foregoing.

Pursuant to 21 CFR 807.39, "Registration of a device establishment or assignment of a registration number does not in any way denote approval of the establishment or its products. Any representation that creates an impression of official approval because of registration or possession of a registration number is misleading and constitutes misbranding."

The U.S. Food and Drug Administration does not issue a certificate of registration, nor does the U.S. Food and Drug Administration recognize a certificate of registration, UCL-REG SERVICE INC. is not affiliated with the U.S. Food and Drug Administration.




UCL-REGSERVICE INC.
602 ROCKWOOD ROAD,WILMINGTON,
NEW CASTLE DE 19802 USA

For and on behalf of
UCL-REG SERVICE INC.

Cert. No.: M20690
Issued Date: 27 March 2020
Expiration Date: 31 December 2020

PRODUCT DETAILS

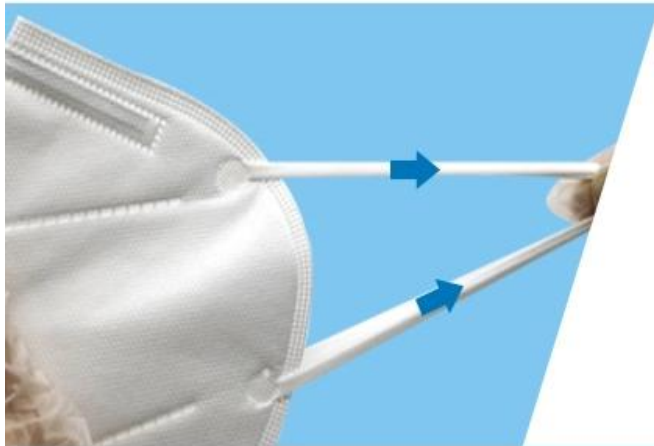

1

High elastic ear loops are comfortable and not tight

2

Adjustable nosepiece helps surgical mask covers the user's face fittingly




3

Solid ultrasonic spot welding edge pressing for strong ear loop and beautiful shape





WEARING A METHOD

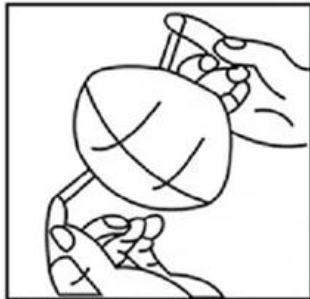

1. Hold the mask so that the stiff bendable strip is on top



2. Pull the bottom of the face mask over your chin.



3. Hold the mask by the ear loops, and place a loop around each ear.



4. Pinch the nose strip to conform to the shape of your nose.



5. Check mask tightness before entering work area .

- Cover your face with both hands and exhale quickly
- If air leaks from the nose, go back to step 4
- If air leaks from the edges, adjust the earband